
AbleNet[®]

iTalk4[™]

www.ablenetinc.com
(800) 322-0956

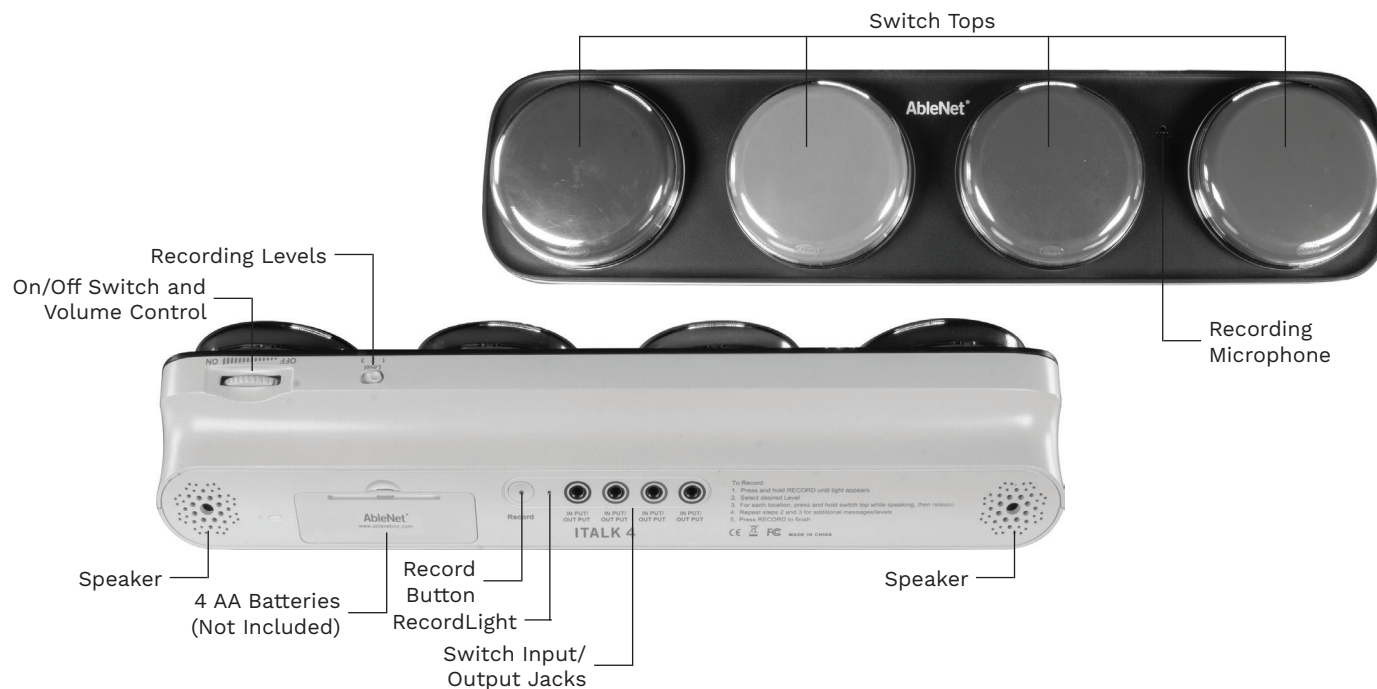Getting Started

1. Insert 4 AA batteries (not included)
2. Turn iTalk4 On
3. Press and hold the record button until the record light activates
4. Press and hold a switch top. Begin speaking into recording microphone when you hear a beep. Release when finished recording message.
5. Repeat for remaining switch tops
6. Press and release record button to exit record mode.
7. Optional: Select a different recording level to record a second and third set of messages.
8. iTalk4 is now ready to use. Press and release desired location to playback message

Recording Levels

iTalk4 provides three recording levels for increased flexibility.

Use the toggle switch on the top of iTalk4 to select the appropriate recording level. You can record separate messages for each message location on each level.

Controlling Toys and Appliances

The iTalk4 switch input/output jacks allow you to activate switch adapted toys or appliance with this device.

1. Plug any switch-adapted toy or appliance into any of the switch input/output jacks
2. Record a message into the corresponding location
3. Press the switch top to play the message and activate the attached toy or appliance. The toy or appliance will turn on for the length of the message.

Mounting

Use the AbleNet Universal Mounting Plate (sold separately) hole pattern D to mount iTalk4 to any of AbleNet's mounting systems. Threaded mounting holes are located on the bottom of iTalk4.

Playing A Message With An External Switch

iTalk4's switch input/output jacks default to controlling toys and appliances. To use external switch(es) with iTalk4:

1. Turn iTalk4 off
2. Attach switch(es) of your choice with 1/8-inch (3.5mm) plugs into the appropriate switch input/output jack(s)
3. While holding any external switch down, turn iTalk4 on

The toy/appliance control feature is disabled when using an external switch(es). To revert to default toy/appliance function, turn iTalk4 off and back on.

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AbleNet, Inc. 2625 Patton Road, Roseville, MN 55113